# MYOMED 632X





### **MYOMED 632X = MYOFEEDBACK**

Myofeedback involves registering the electrical muscle activity using surface electrodes or cavity electrodes. The Myomed 632X is easy to operate and ensures that the therapist and patient receive feedback on the results. The feedback is provided visually via the LCD-screen that clearly displays the pre-set threshold values and the muscle contractions of the patient. One may also opt for an acoustic signal that indicates whether or not the patient achieves the selected values. The rough EMG-signal is audible via headphones. The therapist does not require any labour-intensive supervision, as the Myomed 632X safely guides the patient through the exercises!

### **MYOMED 632X = ALSO PRESSURE FEEDBACK**

In the case of pressure feedback, the Myomed 632X measures changes in pressure by means of cavity electrodes. These changes in pressure come about as a result of contractions of the pelvic muscles. The therapist and patient receive feedback via the LCD-screen.

# **MYOMED 632X = ALSO ELECTROTHERAPY**


The device offers 2 forms of therapy: electro-stimulation and EMG-feedback. These can be used independently or combined. There are various possibilities regarding the combination of EMG-feedback with electro-stimulation: alternate or triggered by a pre-set threshold value.

The Myomed 632X enables the patient to exercise safely, which efficiently promotes recovery!

### **V = ALSO VACUUM**

The vacuum unit (option) is a suction apparatus that is as quiet as a mouse and that offers the possibility of making use of four vacuum electrodes. The suction power can be adjusted across a broad range and can be set both continuous as well as pulsating. Due to a control system, the pump will only be in operation if the vacuum has decreased to below a certain level.

The vacuum electrodes can be positioned quickly and easily. They adapt to the body surface fully and comfortably. Moreover, you can be sure of adequate contact between the body and the entire surface of the electrode.

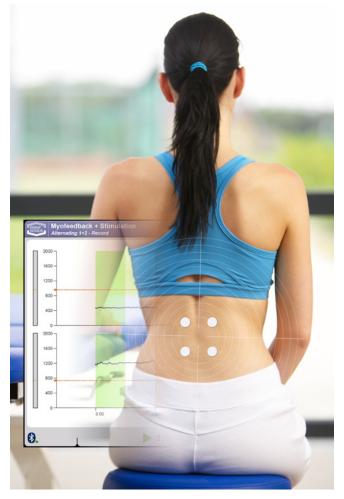

# MYOMED 632X



# 1600966 MYOMED 632 X (1)

Complete device for EMG- and pressure feedback and for electrostimulation



## **1600971 MYOMED 632VX** (1)

 $Complete\ device\ for\ EMG-\ and\ pressure\ feedback\ and\ for\ electrostimulation$ (such as the 1600966), but also fitted with an integrated vacuum unit for easy application of vacuum electrodes for combination and electrotherapy.



# **1600830 MYOCOMBOX**





**ENRAF-NONIUS B.V.** Vareseweg 127 3047 AT Rotterdam The Netherlands
T.: +31-(0)10 - 203 06 00
E: info@enraf-nonius.nl
www.enraf-nonius.com





## X = ALSO COMMUNICATION

The Myomed 632X allows for wireless transmission of data with the Myocombox (optional). The Myocombox can be combined with any random LCD-screen. A larger screen is often desirable in order to improve the visualisation of the myofeedback. The Myocombox makes it possible to display the measurements on a larger screen. Comfortable, particularly for the partially sighted patient, who can view the results without difficulty. The size of the screen is up to you! As the communication is realised by bluetooth wirelessly, it does not matter where the screen and the Myocombox are positioned.


The Myocombox is also used for data transmission. The results of the measurements are sent to the Myocombox, which then programmes the results on an USB-stick. The data can then be processed by an external PC, in a professional report, for example, or in a statement of expenses programme. The Myocombox will prove to be a welcome assistant for any modern therapist!

- Clear feedback displayed on a large computer screen<sup>(1)</sup>
- Wireless data transfer
- Option to export data for evaluation in other software
- Option to create reports

(1) = The Myocombox package does not include the computer screen

# MYOMED 632X, THE MOST COMPLETE DEVICE FOR SAFE TREATMENTS!



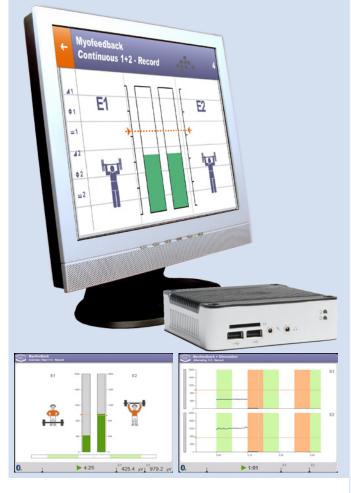

- 2-channel EMG-feedback
- 1-channel pressure-feedback
- 3 Different feedback sessions; continuously, work/rest and template
- EMG- or pressure signal graphically reproduced with adjustable sensibility and time scale
- Continuous reproduction of the EMG or pressure value
- Combined application of EMG-feedback with electrotherapy
- Combined application of EMG- with pressure-feedback
- Special cluster for analysis of the EMG and pressure measurements
- Several surface- and cavity electrodes
- Display of feedback on a large screen possible with the Myocom-
- Wide range of display options
- Data transmission, storage and export possible with the Myocombox

**ENRAF-NONIUS B.V.** Vareseweg 127 3047 AT Rotterdam The Netherlands T.: +31-(0)10 - 203 06 00 E: info@enraf-nonius.nl

# MYOMED 632X



#### STANDARD ACCESSORIES MYOMED (ALL MODELS)

3444673 EMG Cable 6-series

3444674 Cable for reference electrode for EMG-6-series

3444222(1) Set of 8 pcs. adhesive electrodes for stimulation & EMG<sup>(1)</sup> 1460266 Moist pads for rubber electrode 6x8 cm, set of 4

2 x 3444129 Rubber electrodes 6x8 cm, 2 mm female, 2 pc. 2 x 3444211 Patient cable 2-core & 2 mm male plug - black,

with coloured clips

3444020 Strap 100x3 cm Strap 250x3 cm 3444021 3444290 Mains cable 230V-EUR

User manual XXXXXXX

#### **ADDITIONAL STANDARD ACCESSORIES MYOMED V MODELS**

3444078 Dummy plug

2 x 3444507 Vacuum lead hose red 2 x 3444508 Vacuum lead hose black

3444505 Sponges Ø 65 mm, set of 4 (for vacuum

electrodes Ø 60 mm)

2 x 3444503 Vacuum electrodes Ø 60 mm, 2pcs.

(1) = The Myomed is delivered with 1 set of 8 pcs. adhesive electrodes. The articlenumber 3444222 however represents a box of 10 sets.

#### **TECHNICAL SPECIFICATIONS**

EMG-feedback channels: 2, independent

Pressure-feedback:

Combination

EMG/pressure-feedback: yes (2 channel EMG, 1 channel pressure) 0,5-10.000 µV and logarithmically Measurement range EMG:

Measurement range pressure: 0-360 hPa

Current channels: 2, independent, fully galvanically isolated

Current types: 16 current types-36 variants

Programmable positions: unlimited Pre-programmed programmes: >100

Supply voltage: 100-240 V ± 10% (50/60 Hz)

Dimensions (wxdxh): Weight: 6.2 kg

#### OTHER PRODUCTS THAT MIGHT INTEREST YOU

#### **EN-CAR U**

The EN-Car U is your companion for any piece of Enraf-Nonius equipment. For laser, ultrasound, biofeedback, shockwave or electrotherapy, with the EN-Car U you have it all close at hand.

#### **ORDERING DATA**

1468960 EN-Car U

Also see the seperate product leaflets at







# **ACCESSORIES MYOMED 632X**



### **ELECTRODES AND MOIST PADS**

#### **ADHESIVE ELECTRODES**

3444128

| 3444222 | EN-Trode Ø 2,2 cm, 2 mm female, 10 sheets of 8 (also for EMG) |
|---------|---|
| 3444056 | EN-Trode Ø 3,2 cm, 2 mm female,10 sheets of 4                 |
| 3444135 | EN-Trode Ø 5,0 cm, 2 mm female, 10 sheets of 4                |
| 3444057 | EN-Trode 5x5 cm, 2 mm female, 10 sheets of 4                  |
| 3444058 | EN-Trode 5x9 cm, 2 mm female, 10 sheets of 4                  |















| 3444129 | Rubber electrodes 6x8 cm, 2 mm female, set of 2  |
|---------|--|
| 3444130 | Rubber electrodes 8x12 cm, 2 mm female, set of 2 |

Rubber electrodes 4x6 cm, 2 mm female, set of 2.

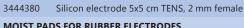

| MOIST TABS FOR ROBBER ELECTROBES |  |
|----------------------------------|--|
| 1460273                          | Moist pads for rubber electrodes 4x6 cm, set of 4  |
| 1460266                          | Moist pads for rubber electrodes 6x8 cm, set of 4  |
| 1460275                          | Moist pads for rubber electrodes 8x12 cm, set of 4 |

| FIXATION STRAPS  |                |
|------------------|----------------|
| 3444020          | Strap 100x3 cm |
| 3444021          | Strap 250x3 cm |
| 3444022          | Strap 100x5 cm |
| 3444023          | Strap 250x5 cm |
| DOINT ELECTRODES |                |



1600820 EN-Point, for making s/d curves

PATIENT CABLES AND ADAPTER PLUGS

|         | Dationt cable 2 care 9: 2 mans made along block |
|---------|---|
| 3444211 | Patient cable 2-core & 2 mm male plug - black,  |
| 3444211 | with coloured clips                             |

### **ADAPTERS**

| 2523524 | Adapter plug, 2 mm female, 4 mm male, red   |
|---------|---|
| 2523523 | Adapter plug, 2 mm female, 4 mm male, black |

### **REMOTE CONTROLS**

| 1497801 | Remote control, 2 channels for 182, 482, 492 and 6-series |
|---------|---|
|         |   |











3444128















1600820











# **ACCESSORIES MYOMED 632X**



#### **MYOFEEDBACK ACCESSORIES EMG ACCESSORIES** 3444673 EMG cable 6-series 3444674 Cable for reference electrode for EMG 3444222 EN-Trode Ø 2,2 cm, 2 mm female, 10 sheets of 8 (also for EMG) 3444117 Anal electrode with fixed cable, 2 mm female plug 3444118 Vaginal electrode 2 mm Evolys, with cable PRESSURE FEEDBACK ACCESSORIES 3444043 Anal pressure sensor CD 3444045 Connecting cable for anal pressure probe 3444043 3444042 Vaginal pressure sensor CD, 6 pcs 3444048 Connection cable for vaginal pressure probe, incl.air inlet $\&\: syringe$ 3444044 Connecting hose for vaginal pressure sensor 3444042 муосомвох











3444118

3444673









3444043





3444044

FOR MYOMED 632VX

Myocombox, for remote lcd display

1600830

| VACUUM A | ACCESSORIES VACOTRON 460                                  |
|----------|---|
| 3444509  | Vacuum electrodes Ø 30 mm, set of 2                       |
| 3444503  | Vacuum electrodes Ø 60 mm, set of 2                       |
| 3444504  | Vacuum electrodes Ø 90 mm, set of 2                       |
| 3444516  | Sponges Ø 30 mm, set of 4 (for vacuum electrodes Ø 30 mm) |
| 3444505  | Sponges Ø 65 mm, set of 4 (for vacuum electrodes Ø 60 mm) |
| 3444506  | Sponges Ø 95 mm, set of 4 (for vacuum electrodes Ø 90 mm) |
| 3444507  | Vacuum lead hose, red                                     |
| 3444508  | Vacuum lead hose, black                                   |
| 3444078  | Dummy plug  |





MYOMD 632X